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Contents

Product introduction	02
User guide	04
Safety precautions	13
Safety precaution include	14
Troubleshooting	16
Technical parameters	17
After service	18

After service, Thank you for choosing the hand trainer.

Please read this manual carefully before using this product.

Please keep this manual properly for guidance when using the rehabilitation training device.

The Hand training device has two modes: the passive mode and Mirror Motion Rehabilitation Training Mode. Our device aims to maintain and expand the range of motion of finger joints, prevent joint contracture and muscle atrophy, greatly promote the reorganization and compensation of the nervous system, and help patients recover as soon as possible.

The recovery of hand dysfunction after stroke is a long-term and arduous process. Patients' conditions, rehabilitation courses, and curative effects are very different. Under normal circumstances, the recovery of hand function is slower than other parts, and long-term rehabilitation training is needed. If you have any questions in use, please get in touch with info@ziglite.com.

If you have any questions in use, please contact: info@ziglite.com

The manual version may upgrade due to changes in software version or technical specifications without prior notice.


Product Introduction

Scope Of Application

It is suitable for rehabilitation training in (families) institutions for patients with finger joint recovery after injury or limb dysfunction sequelae caused by brain injury, stroke, cerebral palsy, spinal cord injury, orthopedic surgery, multiple sclerosis, etc.


Product Composition

Trainer Host

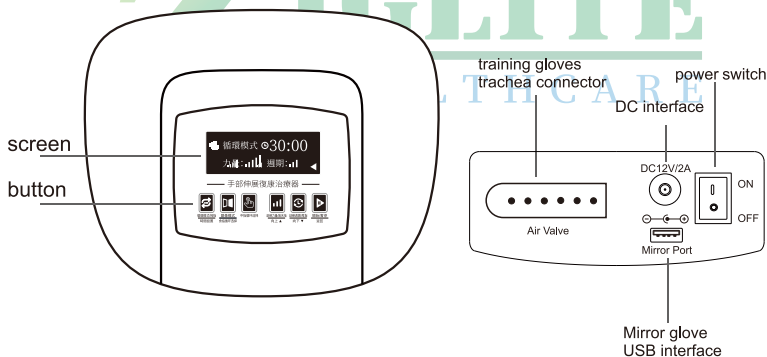
Product Introduction


Training Gloves

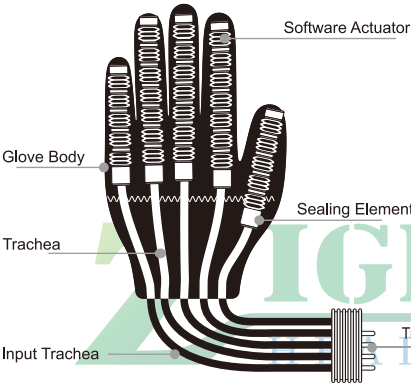Mirror Glove

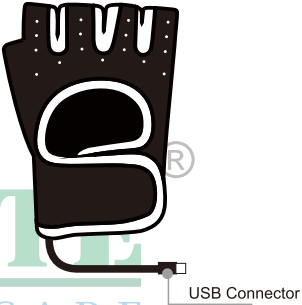
Ower Adapter

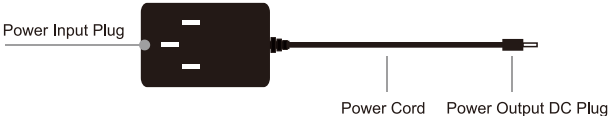User Guide

Operation buttons and interface display

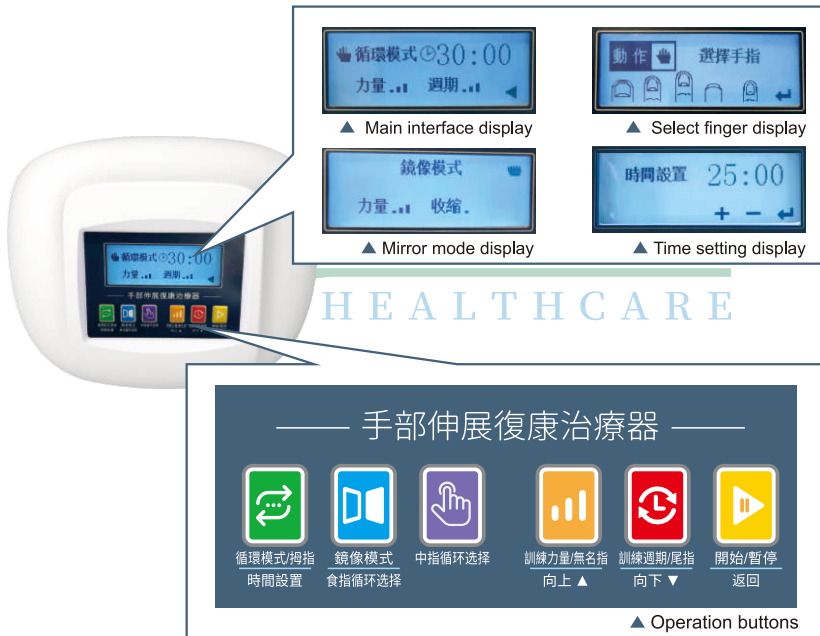User Guide

Detailed button function

**Cycle / split finger mode switching
Thumb training options**

- 1) Press the button once to enter the finger selection interface: corresponding to the finger button to select/stop a finger training.
- 2) In the finger selection interface, press and hold the button for 3s to enter the action mode switch: ①Wave cycle training for each finger/②Cycle training mode switch for all fingers.
- 3) In the main interface, press and hold the button for 3 seconds to enter the training time setting, up to 60 minutes.

**Mirror mode
Index finger training**

- 1) Press the button once to switch the mirroring mode.
- 2) In the finger selection interface, click to select/stop the index finger training exercise.

Middle finger training options

In the finger selection interface, click "Middle finger cycle selection" to select/stop the middle finger training exercise.

**Start/pause
Return**

Start, pause, return and other operations can be performed in the work training mode.

**Cycle time setting
Tail finger training options**

- 1) Press this button to set the cycle time of one grasping and stretching training. The setting value is 1 ~ 5, a total of 5 gears.
- 2) In the finger selection interface, click optional / stop tail finger training.


**Power setting
Ring finger training options**

- 1) Press the "Training Strength" button to set the strength of the training, the set values are: 1, 2, 3, 4, 5, a total of 5 levels of strength.
- 2) In the finger selection interface, click Select/Stop ring finger training.

User Guide

Interface display details

Cycle mode display

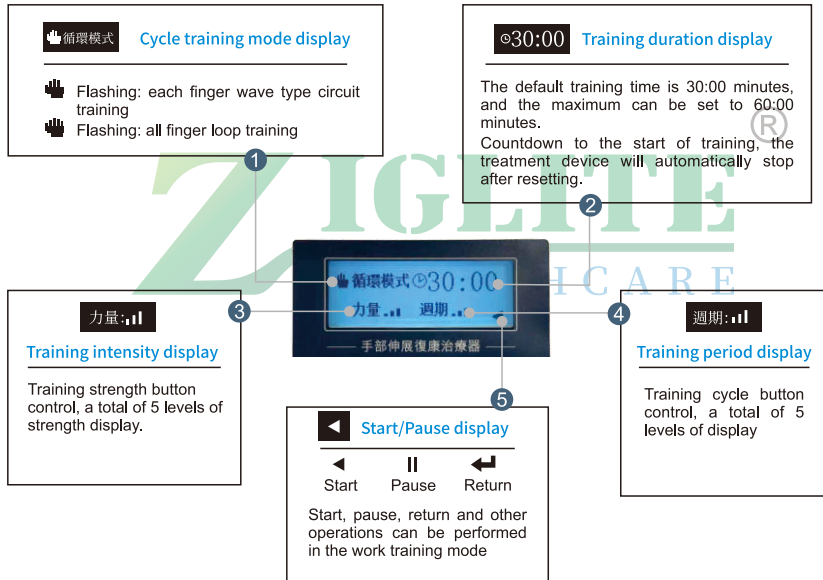User Guide

Interface display details

Mirror mode display

Select finger interface

1 Each finger wave type circuit training mode

Wave circulation mode

Single finger selection / pause (select all finger exercises)

Return to key indication

Flashing means selected, stop flashing means to pause training

2 All finger loop training modes

All finger cycle mode

Select finger training/pause (Choose separate ring finger training)

Return to key indication

Flashing means selected, stop flashing means to pause training

User Guide

Usage of therapy device

1

Place the treatment device host

2

Connect the air pipe to the DC plug

3

Wear training (mirror) gloves

4

Turn on the therapy device

5

Start hand training


6

Therapy is turned off

User Guide

Glove wiring diagram

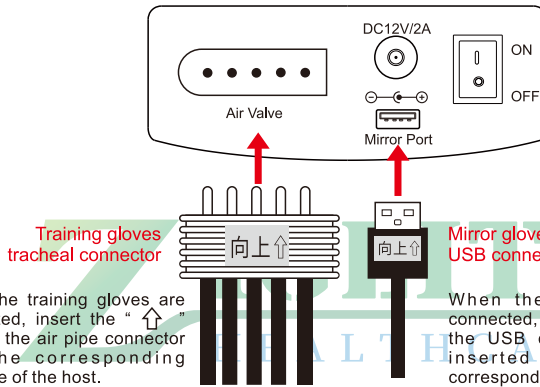When the training gloves are connected, insert the "向上" mark of the air pipe connector into the corresponding interface of the host.

When the mirror glove is connected, the "向上" mark of the USB connector shall be inserted upward into the corresponding interface of the host.

Actual wiring diagram of gloves

User Guide

Cycle training mode

The circulation training mode adopts a fully automatic design. When the patient has no active control ability, it is recommended to use circulation training first, which will help to promote the recovery of flexor and extensor muscle strength of the affected hand, increase the joint's range of motion, and promote the emergence of active finger extension movement.

Default circuit training

	<p>1) Wear training gloves 2) Enter the main training interface after booting</p>		<p>3) After powering on, press the start button to enter the default training: split-finger wave cycle training mode, training time is 30 minutes, cycle strength is the third gear.</p>
--	---	---	--

Set up circuit training as needed

1 Power, period, time setting (4 and 5 can be skipped if you don't need to set the time)

<p>(1)</p> 	<p>(2)</p> 	<p>(3)</p> 	<p>(4)</p> 
<p>Wear training gloves Enter the main training interface after booting.</p>	<p> Set training intensity Press 1 time to set the gear, 1 ~ 5 gears in total.</p>	<p> Set training cycle Press 1 time to set the gear, 1 ~ 5 gears in total.</p>	<p> Set training duration Press and hold for 3 sec to enter the setting,  and  set time plus and minus, each click changes 5 min.</p>
<p>(5) After setting, click  to return to the main interface.</p>			

User Guide

2 Training finger selection (Finger icon flashes for training state, static state is for stop state)

<p>(6)</p> 	
<p>Training finger selection</p> <p>In the main training interface, click  to enter the finger selection interface, select training/stop a finger training.</p>	<p>Finger control: select/stop switching each time you press</p> <p>In this mode, the buttons from left to right correspond to thumb, index finger, middle finger, ring finger and tail finger control.</p>

3 Training loop mode selection

<p>(7)</p>  <p> Finger wave loop</p> <p> Finger circulation</p>	<p>(8)</p>  <p>Return Start training</p>
<p>Training finger selection</p> <p>In the main training interface, long press  for 3 sec to enter the action mode switch:</p> <ol style="list-style-type: none"> Wave type circuit training for each finger Common circulation training mode for all fingers 	<p>Start circuit training</p> <p>After setting, click  to return to the main interface.</p> <p>Click  again to start training.</p>

User Guide

Mirror training mode

The mirror mode refers to the mode in which the patient drives the patient's hand to move synchronously through the exercise of the healthy hand. Through interactive visual feedback, different areas of the patient's brain can be stimulated to produce motor and sensory awareness, and mirror neurons are activated to promote the reconstruction of the patient's motor function and brain function.

<p>(1)</p>  <p>Training glove Mirror glove</p>	<p>(2)</p> 	<p>(3)</p> 	<p>(4)</p> 
<p>Put on gloves and turn on</p> <p>Wear training gloves on the affected hands, Wear mirror gloves for healthy hands, After startup, enter the main training interface.</p>	<p>Enter mirror mode</p> <p>In the main training interface, press the  button once to enter the mirror training mode.</p>	<p>Set training intensity Set stretch duration</p> <p>Press  once each time to set, and the training intensity is 1 ~ 5 gears in total.</p> <p>Press  each time to set the stretching duration, 1-5 gears in total.</p>	<p>Start mirror training</p> <p>The healthy hand grasps (extends) the movement to issue the instruction, and the affected hand grasps (extends) the movement according to the instruction.</p> <p>Note: you need to grasp and press (release) the control button in the mirror glove to start the grip (release) action of the training glove.</p>
	<p>Control the corresponding finger training:</p> <p>There are 5 control keys in the mirror glove, and each key corresponds to 5 fingers of the training glove. The healthy hand sends movement instructions to the corresponding fingers of the training gloves by controlling the buttons in the mirror gloves, and the corresponding fingers of the training gloves receive instructions to make grasping (extension) movement</p>		

Safety Precautions

[warning] failure to comply with the following safety instructions may lead to serious accidents such as fire, electric shock and casualties.

- Do not disassemble, repair, or modify the components of the equipment by yourself to avoid abnormalities in the treatment device.
- Do not store components in or near water heaters, microwave ovens, cooking equipment, or high-pressure containers to avoid fire caused by overheating.
- It cannot be sterilized or autoclaved.
- Do not continue to use it when the alarm prompt appears.
- It is suitable for hand circulation and mirror movement rehabilitation training. Please do not use it for other purposes to avoid accidents.
- Do not store the parts in a place that is too hot or too cold to avoid accidents.
- The validity period of the product is 5 years, and the disposal of the product at the end of its service life should comply with the requirements of local laws and regulations.
- Please keep the power connector dry and water-free, please do not plug or unplug the power supply with wet hands to avoid electric shock or injury.

Working Environment

- 1) Ambient temperature: 10°C~35°C
- 2) Relative humidity: 30%~75%, excluding condensation
- 3) Atmospheric pressure range: 700hpa~1060hpa
- 4) Power supply: AC100~240V, 50/60Hz (main unit power supply of the treatment device DC12V)

Storage And Transportation

- During transportation, please fix the treatment device in a relatively stable place to prevent bumping and damage of the treatment device.
- Please prevent liquid or foreign objects from entering the treatment device to avoid damaging the internal electronic components and causing malfunctions.
- Do not store in a sterile gas or active gas environment to avoid damage to internal electronic components and cause malfunctions.
- Do not fold the gloves forcefully to avoid damaging the soft-actuated actuator and causing malfunctions.
- Do not store the therapy device in a humid environment or strong magnetic field interference and other environments that are not suitable for electronic equipment.

Safety Precaution Include

- In case of skin allergy and other discomfort during the use of the product, please stop using and consult a doctor.
- For patients with diseases or wounds or being treated for them, or with previous medical history, please consult a doctor after being familiar with the following contents.
- In case of sprain, trauma, bone injury of unknown condition, severe osteoporosis, untreated fracture or people who need bracing after treatment, do not use this product without professional confirmation.

Contraindication

- Patients with severe blood flow disorder such as thoracic venous tumor and hemocolumn disease;
- Patients with ulcer, wound, purulent fruit, bleeding, scald, inflammation, etc. on the affected hand. ®

Other Precautions

- 1) In case of abnormal experience or body feeling during use, stop using it and use it under the guidance of a doctor after seeing a doctor;
- 2) In case of redness, pruritus and rash during use, stop use and seek medical attention;
- 3) Do not use it together with other equipment or in the state of applying drugs;
- 4) Do not use it alone for unconscious people or young children;
- 5) Do not use this product for any purpose other than the purpose of use.

Before Use

- 1) Before wearing training gloves, please try to run the therapeutic device host for 1 minute to confirm that all parts function normally;
- 2) Patients with mandatory systolic spasm need to massage before using the therapeutic device;
- 3) Note: before training, please confirm that the non muscle tension of the hand to be trained is too high (above grade II).

Safety Precaution Include

In Use

In case of the following conditions during operation, please stop using immediately and disconnect the power supply.

- 1) The main machine of the trainer makes abnormal sound
- 2) The surface temperature of the trainer is too high
- 3) The main machine of the trainer is flooded
- 4) The user feels pain, swelling or other discomfort symptoms
- 5) If spasm occurs during training, stop using it immediately

Usage Time

- 1) It is suggested to select appropriate training methods according to their own actual situation
- 2) The time of each use should not be too long. If the hand feels sore the next day after training, it is necessary to reduce the number and time of daily training
- 3) It is recommended to use 2-4 times a day for about 20 minutes each time, with an interval of at least 45 minutes
- 4) If you need to change or extend the use time, please consult professional medical personnel

Use Of Power Adapter

- 1) The power adapter must be plugged in and out by the patient's supervisor
- 2) After using the therapeutic device, please cut off the power adapter and 220V power supply
- 3) Do not use a power adapter that is not provided with this machine

Troubleshooting

Problem	The Reason	Measures
Flip the "power switch" down to turn on the machine	Poor power connection	Check the power cord and power adapter, the connection between the power adapter and the host
Training gloves are not working	Gloves that are not properly attached can damage the host or glove structure	Check the physical status of the training gloves, the connection between the treatment device and the training gloves
The applicator gloves cannot be inserted into the applicator host	Incorrect direction of air pipe connector of Training Gloves	Match the arc surface and straight surface of the training glove tracheal joint to the host

Note: If you still can't use the above method, please contact us. 

Care And Maintenance

Clean And Disinfect

- Do not wash training gloves and mirrored gloves directly with water. Dip a soft cloth in water or neutral detergent, wring it out, and wipe the gloves. If necessary, use a rag with a concentration of 75% medical alcohol for cleaning or disinfection, and then use it after it is completely dried.
- When there is dust or dirt on the surface of the main body of the treatment device, the user can wipe it with 75% medical alcohol. When wiping, please use a dry cloth dipped in a small amount of alcohol to avoid dripping or flowing into the product.
- Clean with clean and dry soft cloth, and dry in a cool and ventilated place.

Note

1. The use of cleaning and disinfecting agents other than alcohol may damage the product, shorten the service life or cause safety risks.
2. It is forbidden to use high-temperature drying and exposure to the sun for quick drying or disinfection.
3. Do not use gasoline, thinner or similar organic solvents to clean any part of this product. Do not use radiation, steam, ethylene oxide and other methods for disinfection, otherwise the product will be damaged.

Technical Parameters

Product Specifications

Product Name	Hand Stretching Rehabilitation Therapy Device	Waterproof Level	IPX0
Product Model	SBA-2	Velocity Setting	1-5
Security Type	I Type	Time Setting	5-60 minutes
Power Supply	AC100-240V 50/60Hz	Running Noise	≤60dB
Host Size	251*204*111mm	Training Mode	Cycle training/ Mirror training


Glove Size

Size	S	M	L	XL	XXL
Hand length (mm)	140-160	160-180	165-190	190-220	220-230
Hand width (mm)	60-80	70-90	80-100	100-120	120-130
Middle finger circumference (mm)	45-50	50-60	60-70	70-90	90-100

* The corresponding model can be selected according to the user's actual hand measurement situation.

If two or more measured values are at the critical point, it is recommended to choose one yard larger.

* Women generally choose s/m, men generally choose l/xl, and the larger hand is XXL.

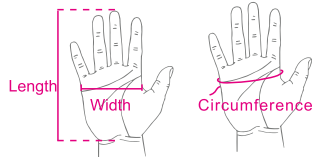After service

After sales email : info@ziglite.com

After sales address : Room B3, 18/F Bonsun Industrial Building,
364-366 Sha Tsui Road, Tsuen Wan, Hong Kong

Packing List		
Number	Name	Quantity
1	Therapeutic device host	1 set
2	Training gloves	1 set
3	Mirror gloves	1 set
4	Power adapter, power cord	1 set
5	Product specification	1 set

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Office Address : Room B3, 18/F Bonsun Industrial Building, 364-366 Sha Tsui Road, Tsuen
Wan, Hong Kong

Email : info@ziglite.com

Fax: +852 3544 0462

Office : +852 3974 5593

Whatsapp : +852 6756 7703

Wechat ID : Gg971213, helen-3110

NO. B02-CN-EN

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Ziglite Smart Health Care Products Co., Ltd.